


EZCare™ Softouch

TRACHEOSTOMY TUBE HOLDER

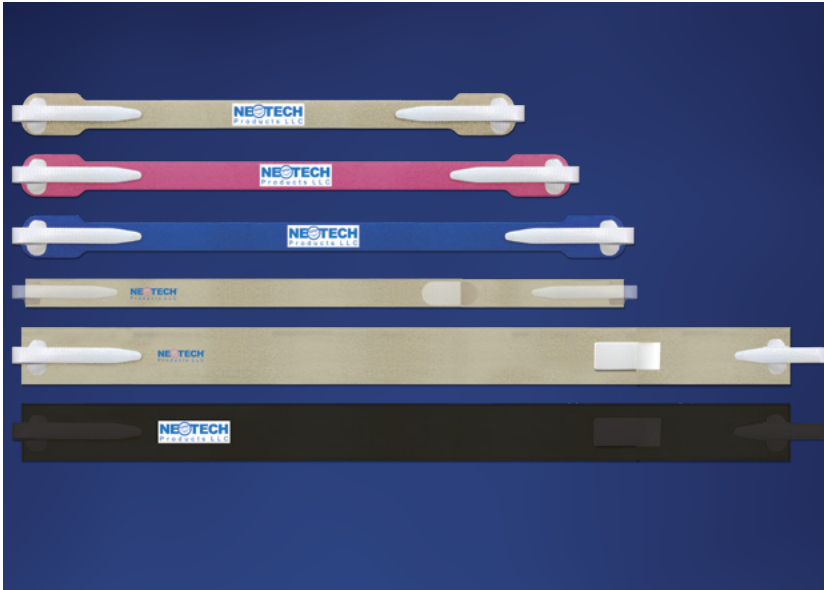

FEATURES & BENEFITS

- Thick NeoFoam® Ultra for better stability
- Soft, breathable, skin friendly NeoFoam wicks perspiration
- Designed for comfort
- Lightweight and form fitting
- N925 and N926 can be cut to size
- Unique material helps to prevent sliding
- Not made with natural rubber latex or plasticizer DEHP
- Individually packaged

CAT. #	SIZE	HCPCS	QTY/
N9106	6" (152mm) Single piece	A7526	20/box
N9107	7" (178mm) Single piece	A7526	20/box
N9108	8" (203mm) Single piece	A7526	20/box
N9109	9" (229mm) Single piece	A7526	20/box
N9110	10" (254mm) Single piece	A7526	20/box
N9111	11" (279mm) Single piece	A7526	20/box
N925*	11" (279mm) 2-piece(tan only)	A7526	20/box
N926**	18" (457mm) 2-piece(tan only)	A7526	20/box
N926BK-10**	18" (457mm) 2-piece(black only)		10/pack

Tan=TN, Pink=PK, Blue=BL

*1/2" (13mm) tall

**1" (25mm) tall



Parent Invented



Made in USA

EZCare™ Softouch

TRACHEOSTOMY TUBE HOLDER

INDICATIONS FOR USE

The EZCare Softouch is intended to secure a patient's tracheostomy tube. It is intended for use on pediatric (neonates, infants and children) patients.

DIRECTIONS FOR USE

Measuring

Select proper size by using EZCare Softouch measuring tape strip provided.

Place start line at flange on left side and wrap behind neck to right flange. If measurement falls between sizes, try larger size first.

Single Piece (N9105-N9111)

- 1 With smooth, white side facing neck, insert hook tab through flange on one side of plate.
- 2 Pull tab back over flange and firmly press white hook fabric to neckband fabric to ensure proper securement.
- 3 Wrap EZCare around neck and insert hook tab through flange on other side of plate.
- 4 Pull tab back over flange and firmly press white hook fabric to neckband fabric to ensure proper securement.

Note: Ensure that EZCare is taught, but not overly stretched. Adjust as needed.

- 5 Replace EZCare Softouch every 5 days or per hospital protocol, whichever is sooner.
Replace immediately if hook tab is worn or visibly soiled.

2-piece (N925-N926)

- 1 Separate two pieces.
- 2 With smooth, white side facing neck, insert hook tabs through flanges on each side of plate.
- 3 Pull tab back over flange and firmly press white hook fabric to neckband fabric to ensure proper securement.
- 4 Wrap both sides of EZCare around neck and bring together
- 5 Determine proper length, pull longer piece away from neck and cut excess.

Note: Ensure that skin will not be exposed to hook fabric secure tab.

- 6 Bring both sides back together and firmly press white hook fabric secure tab to neckband fabric to ensure proper securement.

Note: Ensure that EZCare is taught, but not overly stretched. Adjust as needed.

- 7 Replace EZCare Softouch every 5 days or per hospital protocol, whichever is sooner.

Replace immediately if hook tab is worn or visibly soiled.

GENERAL INFORMATION

Product Materials

NeoFoam Ultra, hook fabric

This product is not made with natural rubber latex or plasticizer DEHP.

Size

N9105	5" (127mm)	Single piece
N9106	6" (152mm)	Single piece
N9107	7" (178mm)	Single piece
N9108	8" (203mm)	Single piece
N9109	9" (229mm)	Single piece
N9110	10" (254mm)	Single piece
N9111	11" (279mm)	Single piece
N925	11" (279mm)	2-piece (tan only)
N926	18" (457mm)	2-piece (tan only)
N926BK-10	18" (457mm)	2-piece (black only)

Packaging

Individually packaged: Yes

Package type: Polyethylene bag

Dispenser box: Yes

Visit neotechproducts.com for the most up-to-date information.

Boxes are printed with environmentally friendly ink on recyclable material.

CAUTIONS & WARNINGS

Federal law restricts this device to sale by or on the order of a physician.

Replace immediately if hook tab is worn or visibly soiled. Worn tab may cause decannulation.

Check regularly to prevent disengagement.

Discontinue immediately if skin irritation occurs.

Single patient use only.

Certified



Corporation

For more information or to order your FREE sample call
1-800-966-0500 or visit neotechproducts.com

Neotech Products LLC • 28430 Witherspoon Pkwy, Valencia, CA 91355



Made in USA